

BioArk® 100bp DNA Ladder

Cat. # : BADM3367

Product Information

Product Name	Cat. No.	Spec.
GN100bp DNA Ladder III	G3367-01	500 µL
	G3367-05	5×500 µL

Product Description/Introduction

The GN100bp DNA Ladder III consists of 10 linear double-stranded DNA (100 bp, 200 bp, 300 bp, 400 bp, 500 bp, 600 bp, 700 bp, 800 bp, 900 bp, 1000 bp) bands mixed with 1 x DNA Loading Buffer containing a blue indicator dye and is suitable for the analysis of DNA bands in agarose gel electrophoresis. This product is ready-to-use, according to the experimental needs, take 5-10 µL electrophoresis, easy to use, clear electrophoresis bands, easy to accurately determine the content of the target product DNA. The 500 bp band in GN100bp DNA Ladder III has a DNA concentration of 100 ng/5 µL and shows a bright band; the remaining bands have a DNA concentration of 50 ng/5 µL.

Storage and Shipping Conditions

Ship with wet ice; store at -20°C, valid for 24 months.

Product Contents

Component	G3367-01	G3367-05
GN100bp DNA Ladder III	500 µL	5×500 µL
10 x DNA Loading Buffer (Peach)	1 mL	1 mL
Manual	One copy	


Assay Protocol / Procedures

1. After thorough thawing and mixing, take 5 µL of the product and add it to the sampling wells of the agarose gel (the volume of the sample volume is proportional to the width of the wells: 1 µL: 1 mm, if the sampling wells are wider, the volume of the sample volume can be increased appropriately), and then carry out electrophoresis.
2. Recommended electrophoretic conditions are 2% agarose (G5056 recommended) gel at 4-10 V/cm.
3. Stain by EB or other nucleic acid dyes (G3606 recommended) and observe the electrophoretic bands under a UV lamp or gel imager.

Note

1. When using, thaw and mix thoroughly, thaw and store at 4°C, avoid freeze - thaw cycles.
2. Recommended electrophoresis buffer: TAE/TBE (recommended G3001, G3002); change the electrophoresis buffer and use freshly prepared agarose gel in time to avoid affecting the electrophoresis results.
3. For your safety and health, please wear safety glasses, gloves, or protective clothing.

GN100bp DNA Ladder III

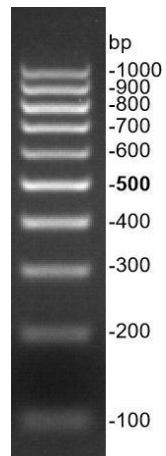

2 % TAE agarose gel stained with EB (loading with 5 μ L DNA marker)

For Research Use Only!