

BioArk® Prestained Protein Marker 250KD

Cat. #: BAPM1126

Product Information

Product Name	Cat.No.	Spec.
Prestained Protein Marker 250KD	BAPM1126-01	1 mL
	BAPM1126-02	5×1 mL
	BAPM1126-05	15×1 mL

Product Description:

Western Protein Marker 250KD consists of eight high-purity recombinant proteins and two pre-stained recombinant proteins, which is convenient for monitoring of electrophoresis and transfer efficiency. It is suitable as a protein molecular weight standard for SDS-PAGE and Western Blot.

The bands display in Tris-Glycine gel range from 10-250 kDa (~10, ~15, ~25, ~35, ~55, ~70, ~100, ~130, ~250 kDa), with orange-red band at 70 and 25 kDa, green band at 10 kDa, and the others can bind almost all types of antibodies (except chicken antibodies). Western Protein Marker 250KD can bind the antibody to the target protein at the same time and develop color by ECL or other methods.


Key Features:

(~10, ~15, ~25, ~35, ~55, ~70, ~100, ~130, ~250 kDa)
orange-red band at 70 and 25 kDa, green band at 10 kDa

Delivery and storage:

Ship with wet ice; Store at -20°C; Valid for 12 months.

BAPM1126

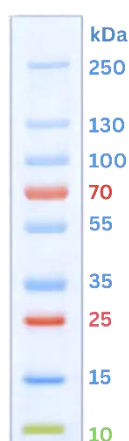

Prestained Protein Marker 250KD

For Research Use Only!